# concept

Osobní váha  
Osobná váha  
Waga osobowa

Personiskie svāri  
Personal scale  
Personenwaage

**VO2850**

**CZ**

**SK**

**PL**

**LV**

**EN**

**DE**

## PODĚKOVÁNÍ

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## ACKNOWLEDGMENT

Thank you for purchasing a Concept product. We hope you will be satisfied with our product throughout its service life.

Please study the entire Operating Manual carefully before you start using the product. Keep the manual in a safe place for future reference. Make sure other people using the product are familiar with these instructions.

| TECHNICAL PARAMETERS |                                                                 |
|----------------------|-----------------------------------------------------------------|
| Voltage              | DC 9 V, power supply 1 piece battery 6LR61 (alkaline, supplied) |
| Weighing range       | 5 - 160 kg                                                      |
| Weighing increments  | 0,1 kg                                                          |
| Weighing units       | kg, st, lb                                                      |
| Age setting range    | 10 - 99 years                                                   |
| Height setting range | 100 – 220 cm                                                    |

### WARNING!!!

**The unit is not for use by pregnant women or people with pacemakers!**

**When measuring the body fat and water of persons with metal implants or in those who are taking contraceptives, the results can be distorted or erroneous.**

## IMPORTANT SAFETY INSTRUCTIONS:


- Use the unit only as described in this operating manual.
- Remove all the covering and marketing materials from the unit before the first use.
- Make sure you follow the marked polarity of the batteries when inserting and replacing them.
- If you are not going to use the unit for an extended period of time, remove the batteries.
- Do not leave the unit unattended when in use.
- Do not allow children or unskilled people to handle the unit. Use the unit out of the reach of these individuals
- Persons with limited movement capacity, reduced sensory perception, insufficient mental capacity or those who are unaware of the proper handling should use the product only under the supervision of a responsible person aware of these instructions.
- Pay the utmost attention when using the unit close to children.
- Don't let anyone use the appliance as a toy.
- Use the unit only on a stable, horizontal surface.
- Before cleaning, assembly, or disassembly of the cover, switch the unit off and remove the batteries!
- Do not use abrasive or chemically aggressive substances to clean the unit.
- If segments or digits are incorrectly displayed, try the so-called reset. Remove the battery, wait for a couple of seconds and put it back. If this procedure does not help, contact the vendor or an authorised service centre.
- Do not use the unit if it does not operate properly, or if it has been dropped, damaged or exposed to any liquids. Take the unit to an authorised service centre for testing and repair.
- The appliance is intended for domestic use only, not for commercial use.
- Keep the appliance away from any other heat sources such as radiators, ovens etc. Do not expose the unit to direct sunlight or moisture.
- Do not touch the unit with wet or damp hands.

- The product is a fragile measuring device. Prevent it from crashing or falling. Do not jump on the weighing surface.
- Do not repair the appliance yourself. Refer to an authorised service centre instead.


**Failure to follow the manufacturer's instructions may lead to refusal of warranty repair.**

## PRODUCT DESCRIPTION

1. Weighing surface
2. Display
3. MODE button
4. Buttons ▼ and ▲
5. Measuring electrodes
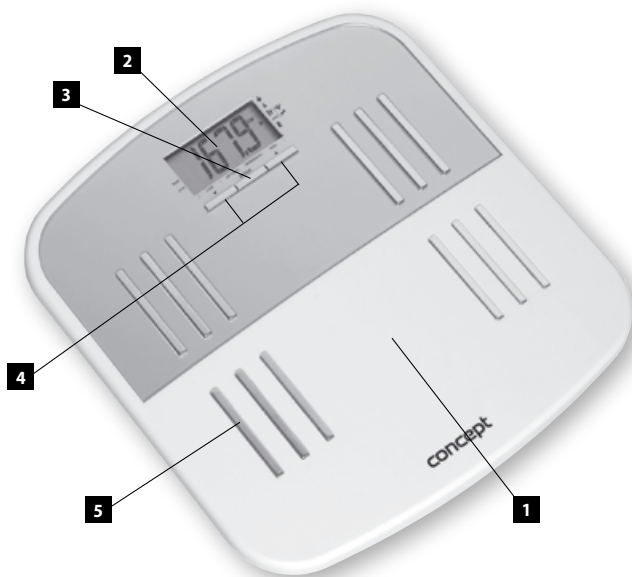
## DESCRIPTION OF THE DISPLAY
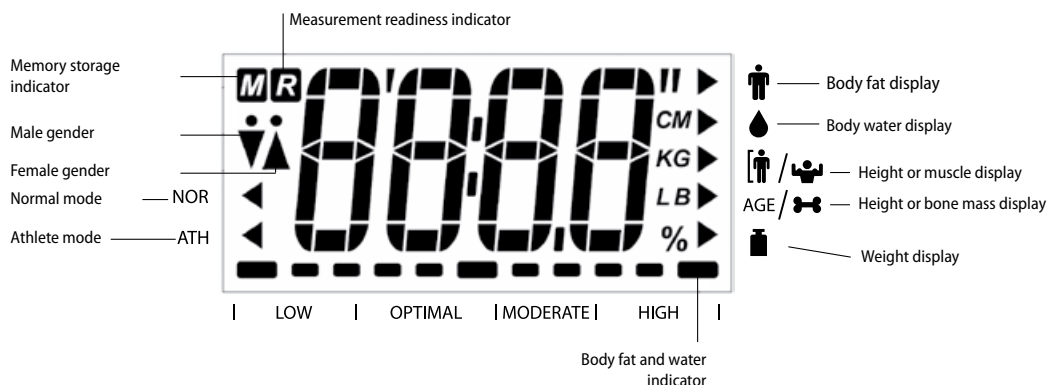
## FUNCTION OF MEASURING BODY FAT AND WATER

The unit works on the principle of measuring bioelectric resistance using a weak electric current from measuring electrodes located on the weighing surface. This electric current is absolutely harmless; the user does not feel it flowing from one leg to the other.

The water contained in non-fat body mass conducts better than the fat layer. The total measured resistance, bioelectric impedance, allows calculation of the percentage of fat and non-fat tissue with regard to weight, age, gender and height.

### RECOMMENDATION

Weigh yourself on a stable horizontal surface. Use a hard surface, not a carpet.

Weigh yourself on the same spot. An uneven surface may affect the results.

Weigh yourself without clothes, before a meal, at approximately the same time of day; ideally, in the morning.

When measuring body fat and water you have to be barefoot, and the soles of your feet must be dry.

## OPERATING INSTRUCTIONS

### INSERTING THE BATTERIES

1. Open the battery compartment cover located underneath the unit (Fig. 1).
2. Insert one battery type 6 LR61 (9V alkaline, supplied) in the battery compartment, and replace the cover carefully.
3. after every battery replacement calibration of the unit is necessary: Put one foot in the middle of the measuring surface to make the unit switch on. The „0000“ message will be displayed (Fig. 2). Wait until the unit switches off.

### SELECTION OF WEIGHING UNITS

When switched off, set the required units (lb, st-lb, kg) using the switch located on the underside (Fig. 3).

### WEIGHING

1. Put the unit on a stable, horizontal surface.
2. Step onto the weighing surface carefully (Fig. 4). Make sure not to slide, or turn the unit over. The result will appear on the display (Fig. 5).
3. Step off the weighing surface.
4. The unit will switch off automatically after a couple of seconds.

### SETTING PERSONAL DATA IN THE MEMORY

You may save the personal data of up to 4 people in the unit's memory.

1. Put the unit on a stable, horizontal surface.
2. Push the MODE button.
3. Select your height in cm using the arrows (Fig. 6). Confirm by pushing the MODE button.
4. Select your age using the arrows (Fig. 7) and confirm by pushing the MODE button.
5. Select your gender using the arrows (Fig. 8) and confirm by pushing the MODE button.

Fig. 1

Fig. 2

Fig. 3

Fig. 4

Fig. 5

6. Select the normal or athlete mode using the arrows (Fig. 9) and confirm by pushing the MODE button. A sound signal will be heard and the readiness indicator displayed.

In case you want to store the data in the unit's memory, select the required memory number using the arrows (Fig. 10) and push the MODE button. Occupied memory is indicated by the symbol „■“.

7. A sound signal will be heard and the memory storage indicator (M) and measuring readiness indicator (R) will be displayed.

8. You can start to measure the body fat, water, muscle and bone mass.

### ERASING PERSONAL DATA FROM THE MEMORY

1. Push the MODE button for 2 seconds when the unit is switched off.
2. The memory number and entered personal data will be alternately displayed.
3. Select the required memory number using the arrows.
4. Press both the arrows simultaneously. The selected memory will be erased.

### BODY FAT, WATER, MUSCLE AND BONE MASS MEASURING FUNCTION

1. Push some of the arrows or the MODE button for 2 seconds when the unit is switched off.
2. The memory number and entered personal data will be alternately displayed.
3. Select the required memory number using the arrows and confirm by pushing the MODE button.. A sound signal will be heard and the memory storage indicator (M) and measuring readiness indicator (R) will be displayed. Step on the weighing surface barefoot; place your feet on the measuring electrodes (Fig. 4). Make sure not to slide, or turn the unit over. The resulting weight will be displayed (Fig. 5).
4. After a while the measurement of body fat and water will be displayed (Fig. 11). To ensure a correct measuring the feet, shins and thighs shall not be in contact.
5. After measuring, the value of body fat will be displayed (body weight percentage, Fig. 12), together with the value of water content (body weight percentage, Fig. 13), body muscle mass (body weight percentage, Fig. 14), and the bone mass (body weight percentage, Fig. 15).
6. The unit will switch off automatically after a couple of seconds.

### SWITCHING THE UNIT OFF

You can switch the unit off by pushing the MODE button for 2 seconds.

The unit will switch off automatically if you don't push any buttons for 10 seconds.

Fig. 6

Fig. 7

Fig. 8

Fig. 9

Fig. 10

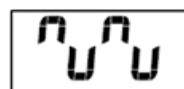
Fig. 11

Fig. 12

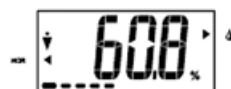
Fig. 13

Fig. 14

Fig. 15

## OTHER INFORMATION ON THE DISPLAY

| Displayed data       | Cause and solution                                                                                                                           |
|----------------------|----------------------------------------------------------------------------------------------------------------------------------------------|
| Nothing is displayed | Weak or flat battery, replace with new                                                                                                       |
| Lo                   |                                                                                                                                              |
| Err0                 | Calibration error. Calibrate the unit again.                                                                                                 |
| Err1                 | Measuring error. Repeat the process, stand still during the measuring.                                                                       |
| Err2                 | The unit has been overloaded and can be damaged.                                                                                             |
| Err3                 | Body impedance measurement error. Make sure you stand on the measuring electrodes correctly, barefoot. If necessary wet your soles slightly. |
| Err4                 | Body fat measurement error. If necessary wet your soles slightly.                                                                            |
| Err5                 | Body water measurement error. If necessary wet your soles slightly.                                                                          |

## INDICATIVE TABLE OF VALUES AND THEIR MEANING

**Note:** The tables and values serve for informational purposes only, so that you may observe long-term changes. The measured value can differ from the actual value, as the accuracy of the unit depends on humidity, the skin of the soles of the feet and other bodily influences (e.g. menses). Contact your physician for more detailed information. Muscles of the same volume are heavier than fat. The weight of clothing is added to the fat weight.

### Body fat percentage

| Females |             |               |            |       |
|---------|-------------|---------------|------------|-------|
| Age     | Underweight | Normal weight | Overweight | Obese |
| 20 - 29 | < 19        | 20 - 28       | 29 - 31    | > 31  |
| 30 - 39 | < 20        | 21 - 29       | 30 - 32    | > 32  |
| 40 - 49 | < 21        | 22 - 30       | 31 - 33    | > 33  |
| 50 - 59 | < 22        | 23 - 31       | 32 - 33    | > 34  |
| 60 +    | < 23        | 24 - 32       | 33 - 35    | > 35  |

### Body fat percentage

| Males   |             |               |            |       |
|---------|-------------|---------------|------------|-------|
| Age     | Underweight | Normal weight | Overweight | Obese |
| 20 - 29 | < 13        | 14 - 20       | 21 - 23    | > 23  |
| 30 - 39 | < 14        | 15 - 21       | 22 - 24    | > 24  |
| 40 - 49 | < 16        | 17 - 23       | 24 - 26    | > 26  |
| 50 - 59 | < 17        | 18 - 24       | 25 - 27    | > 27  |
| 60 +    | < 18        | 19 - 25       | 26 - 28    | > 28  |

## Body water percentage

|       | Woman  | Man    |
|-------|--------|--------|
| Age   | Normal | Normal |
| 6-12  | 58-70  | 60-73  |
| 13-79 | 42-61  | 46-66  |

**Athlete mode** - Athletes have a higher percentage of muscle, and a higher tendency to dehydrate, than normal adults. They would have a higher body fat result in normal mode. The athlete mode is designed for people who do sports intensively, at least 6 hours a week.

**Portion of body muscle mass** - this function shows the weight of both striated and smooth muscles, including the water in them.

**Bone mass** - means the volume of calcium and other minerals in the body. Researches have proven that exercise and development of bone tissue helps to develop stronger and healthier bones. As the bone structure does not change in the short term, it is important to develop and keep healthy bones through healthy food and sufficient activity.

## CLEANING AND MAINTENANCE

The surface of the unit may only be cleaned using a wet cloth.

Do not use any detergents or hard objects as they may damage the unit's surface!


**Never clean the unit under running water, do not rinse it, and do not submerge it in water!**

## REPAIR & MAINTENANCE

Any extensive maintenance or repair that requires access to the inner parts of the product must be performed by an expert service centre.

## ENVIRONMENTAL PROTECTION

- Recycle packaging materials and old appliances rather than throwing them away.
- The appliance box may be disposed of at a sorted waste collection point.
- Deliver polyethylene (PE) plastic bags to a waste collection facility for recycling.
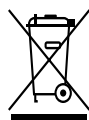
### Recycle the appliance at the end of its service life:

A symbol on the product or its packaging indicates that this product should not be disposed of as household waste. It must be taken to a collection point authorised for recycling electric and electronic equipment. By making sure this product is disposed of properly, you will help prevent negative effects on the environment and human health that would otherwise result from inappropriate disposal of this product. You can learn more about recycling this product from your local authorities, a household waste disposal service, or in the store where you purchased this product.

Disposal of accumulators and battery cells

If an appliance contains an accumulator or a battery cell, these components contain substances that may damage the environment when decomposing. Prior to the appliance disposal or the accumulator or battery cell replacement, it is necessary to remove the original component and take it to a suitable designated collection yard.

### WARNING!

#### An accumulator or battery cell:

- **Never throw into a fire!**
- **Never disassemble it!**
- **Never short-circuit its terminals!**

Should an accumulator or battery cell leakage occur, the leaked solution is of an alkaline character. If your skin gets stained, rinse it immediately with soap and water and seek medical treatment.

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