## **Berger TB Optimal Sonic Toothbrush**

User`s Manual

[EN]

## Safety Precautions

When using your Rechargeable Sonic Toothbrush, basic safety precautions should always be followed, including:

- Do not open, repair, or alter your toothbrush. Contact a qualified service centre if your toothbrush is in need of service.
- Do not drop or crush your toothbrush. If damaged or defective do not use.
- To ensure safety, only use the parts that come included with your toothbrush.
- If your teeth bleed after using the toothbrush or bleeding continues for more than a week, consult your dentist immediately.
- If you have recently had gum or oral surgery, consult your dentist before using this product.
- Parental supervision is advised for children who use this toothbrush.

 <p><b>CAUTION</b></p>	<p><b>BATTERY WARNING:</b></p> <ul style="list-style-type: none"><li>-The device is equipped with an integrated lithium ion battery. The battery cannot be replaced. Do not attempt to remove the battery from the device.</li><li>-Do not dispose of in fire or expose to excessive heat.</li><li>-Do not crush, puncture, incinerate, or short circuit external contacts.</li><li>-Please recycle or dispose of the battery properly. Contact your local recycling facilities and/or the manufacturer for further information.</li></ul>
---	--

## Introduction

Thank you for purchasing the Berger Rechargeable Sonic Toothbrush, model TB Optimal. Your toothbrush is designed not only to ease sensitive teeth but to also polish and remove stains. In order to get the best results, read this manual carefully before using your new toothbrush.

<u><b>Package Contents</b></u>	<u><b>Specifications</b></u>
<ul style="list-style-type: none"><li>- Rechargeable Sonic Toothbrush</li><li>- Protective Cap</li><li>- Two Replacement Brush Heads</li><li>- User's Manual</li></ul>	<p><b>Charging Time:</b> ~8 hours <b>Working Temperature:</b> 0-50°C <b>Voltage:</b> DC 3.7V <b>Charging Power:</b> 1W</p>

## Features

### - USB Rechargeable

*LED Indicator with low battery warning & power-off protection 3+ weeks of brushing power on a Single Charge.*

### - Oscillating Head

*40,000 Strokes per minute effectively removes plaque without irritating your gums.*

### - IPX7 Waterproof



*Perfect for bathroom use; easy to clean under the faucet.*

### - Smart Timer

*Two Minute auto-timer ensures dentist recommended brushing time.*


*30 Second vibration intervals remind you to thoroughly brush each oral quadrant.*

## A Closer Look

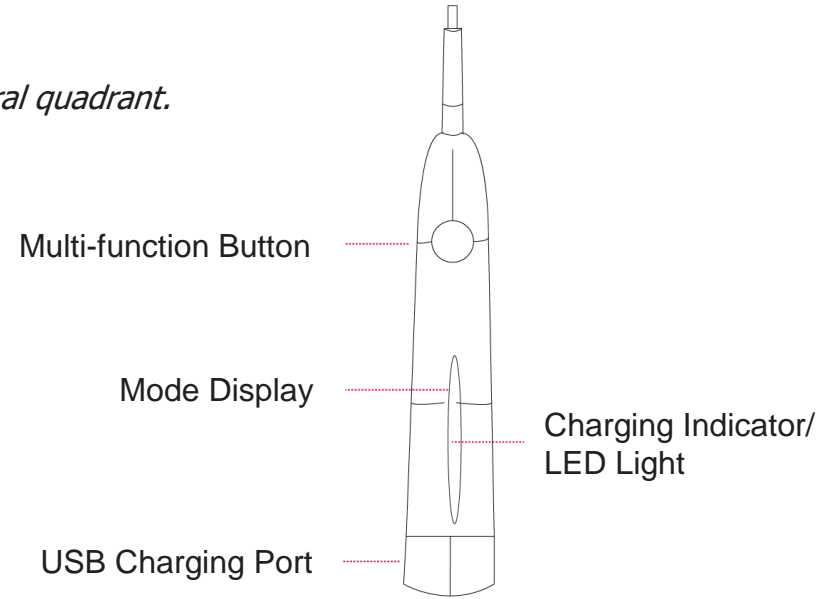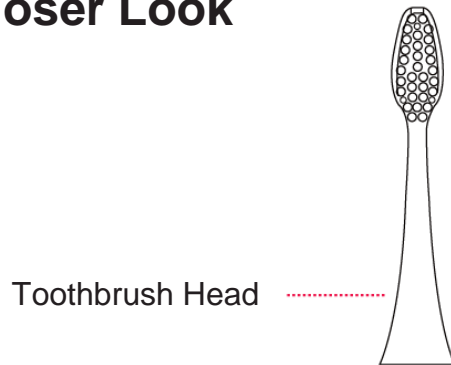
## Attaching the Brush Head

### Install

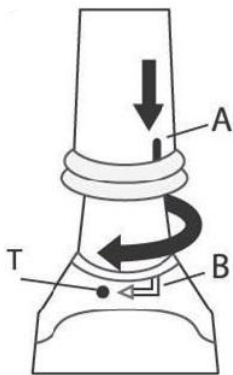Slide the toothbrush head down and turn clockwise so A meets B.

### Detach

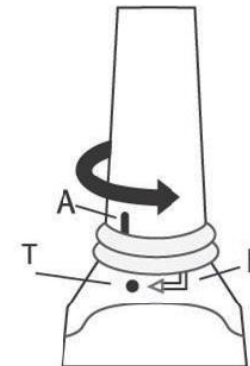Turn the toothbrush head counter clockwise until A meets B. Then pull up to remove.

## Charging the Sonic Toothbrush

1. Connect a Micro USB cable to your toothbrush's charging port.
2. The battery life of your toothbrush will be indicated with the charging indicator light:
  - All 5 LED lights will illuminate with a battery level >80%
  - The bottom 4 white LED lights will illuminate with a battery level at 70-80%
  - The bottom 3 white LED lights will illuminate with a battery level at 60-70%
  - The bottom 2 white LED lights will illuminate with a battery level at 40-60%
  - The bottom white LED light will illuminate with a battery level at 20-40%
  - A red LED light will flash with a battery level at 10-20%
  - Once your toothbrush's battery level falls below 10%, the toothbrush will power off.
3. The indicator light will be on while charging. When fully charged, the indicator light will turn off automatically. **It will take approximately eight hours to fully charge.**
4. If your toothbrush has not been charged for a long period of time, it may take several minutes for the charging lights to illuminate.

## Using the Sonic Toothbrush

1. Turn the toothbrush on. Wet the bristles of the toothbrush and apply the appropriate amount of toothpaste.
2. Select the preferred mode:
  - White:** Recommended to Remove Surface Stains & Plaques
  - Clean:** Recommended for daily use, standard everyday cleaning mode
  - Sensitive:** Recommended for sensitive teeth & gums
  - Polish:** Recommended for use to increase shine on teeth
  - Massage:** Recommended for individuals with sensitive gums
3. Gently move the brush head across your teeth for 10 seconds at a time.
  - It is recommended to go through sections of your mouth when the auto reminder briefly pauses every 30 seconds.
4. When you are finished, hold down the power button for 3 seconds and then immediately let go of the power button to turn the toothbrush off.

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