
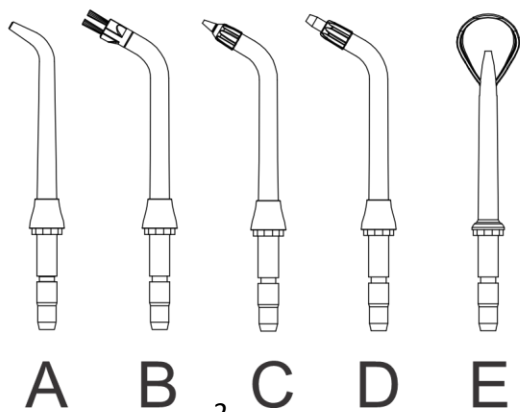
## AD 2185

EN	User Manual	3	PL	Instrukcja obsługi	69
DE	Bedienungsanleitung	6	IT	Manuale d'uso	73
FR	Mode d'emploi	11	SV	Bruksanvisning	78
ES	Manual de usuario	15	BG	Ръководство за употреба	82
PT	Manual de Instruções	19	DA	Brugsanvisning	87
LT	Naudojimo instrukcija	23	SK	Používateľská príručka	91
LV	Lietošanas instrukcija	27	BS	Korisnički priručnik	95
ET	Kasutusjuhend	31	MK	Упатство за употреба	98
HU	Használati útmutató	35	HR	Korisnički priručnik	103
RO	Manual de utilizare	39	UK	Керівництво користувача	107
CS	Návod k použití	43	SR	Упутство за употребу	111
RU	Руководство пользователя	47	AR	دليل المستخدم	115
EL	Εγχειρίδιο χρήσης	52	AZ	İstifadəçi təlimatı	117
NL	Gebruikershandleiding	57	SQ	Manuali i përdorimit	121
SL	Navodila za uporabo	61	KA	ინსტრუქცია	125
FI	Käyttöohje	65			
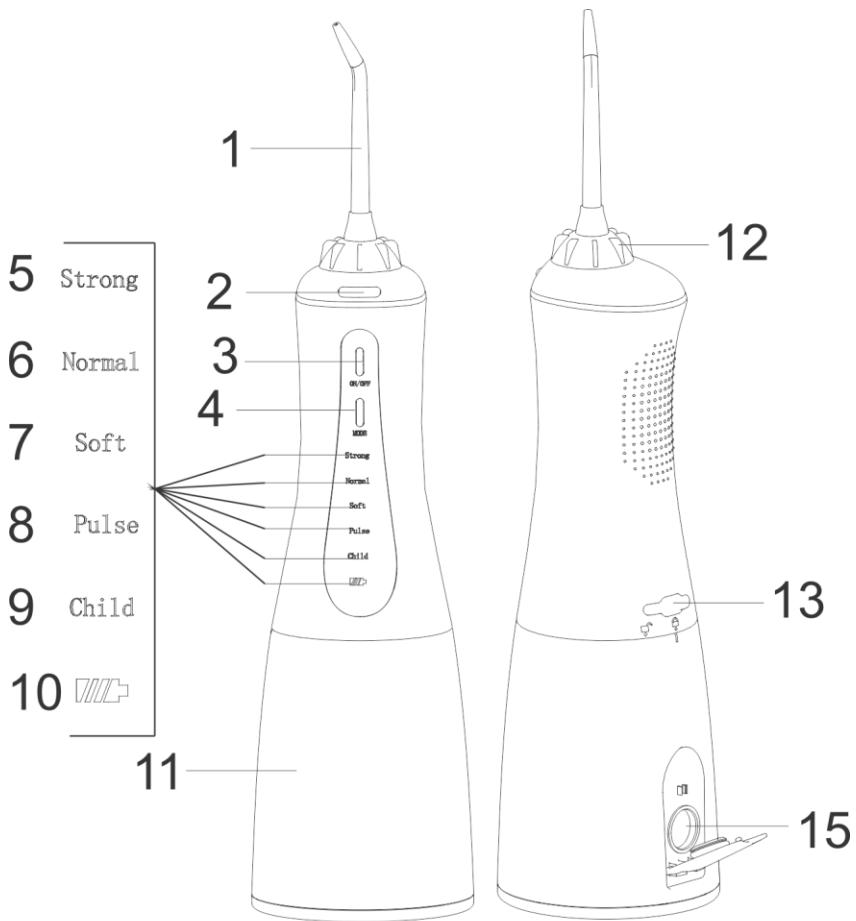
# User Manual (EN)

## GENERAL SAFETY CONDITIONS IMPORTANT SAFETY INSTRUCTIONS FOR USE READ CAREFULLY AND KEEP FOR FUTURE REFERENCE

1. Before using the device, read the instruction manual and follow the instructions contained therein. The manufacturer is not responsible for damages caused by using the device contrary to its intended use or improper operation.

2. The device is for domestic use only. Do not use for other purposes incompatible with its intended use.

3. Connect the device only to a power source of 5V DC 1A or 5V DC 2A .

In order to increase the safety of use, multiple electrical appliances should not be connected to a single power circuit at the same time.

4. Be especially careful when using the device when children are around. Do not allow children to play with the device do not allow children or anyone unfamiliar with the device to use it.

5. **WARNING:** Toothbrushes, may be used by children over 8 years of age and persons with limited physical, sensory or mental ability, or persons without experience or familiarity with the equipment, if it is done under the supervision of a person responsible for their safety or they have been given instructions on the safe use of the equipment and are aware of the dangers of its use. Children should not play with the equipment. Cleaning and maintenance of the equipment should not be performed by children unless they are over 8 years of age and these activities are performed under supervision.

6. **WARNING:** Irrigators, may be used by children over 8 years of age and persons with limited physical, sensory or mental ability, or persons without experience or knowledge of the equipment, if done under the supervision of a person responsible for their safety.

7. **WARNING:** Do not leave the device plugged into an outlet unattended.
8. **WARNING:** Use oral irrigators only with cold or warm water or special solutions, if specified in the manufacturer's instructions.
9. **WARNING:** The device should only be used with the power adapter supplied with the device.
10. **WARNING:** Whenever charging is complete, always remove the plug from the power outlet by holding the outlet with your hand. **DO NOT** pull on the power cord.
11. **WARNING:** Do not use the device near flammable materials.
12. **WARNING:** Do not grip the power adapter with wet hands.
13. **WARNING:** When charging the battery, the equipment must be used so that, at the extreme position of the cord, it is out of reach of a person using a bath or shower.
14. Do not immerse the cable, plug and entire power supply in water or other liquid. Do not expose the device to weather conditions (rain, sun, etc.)
15. Periodically check the condition of the power cord. If the power cord is damaged, the equipment must not be used. A damaged cord should be replaced by a specialized repair facility to avoid danger.
16. Do not use the equipment with a damaged power cord or if it has been dropped or damaged in any other way or is not working properly. Do not repair the device yourself, as there is a risk of electrocution. Return a damaged device to a competent service center for inspection or repair. All repairs may be made only by authorized service centers. Improper repair may cause serious danger to the user.
17. The power cord should not touch heated parts of the device and should not be near other heat sources.
18. To provide additional protection, it is advisable to install in the electrical circuit, a residual current device (RCD) with a rated residual current of not more than 30 mA. In this regard, consult a specialist electrician.

19. Do not use aggressive detergents in the form of emulsions, lotions, pastes, etc. to wash the housing, as they may, among other things, remove the applied informative graphic symbols such as graduations, markings, warning signs, etc.

20. The equipment is not designed to work with external timers or a separate remote control system.

Device description Tooth irrigator:

- |                          |  |                                |
|--------------------------|--|--------------------------------|
| 1. Nozzle                | 2. Nozzle release button                   | 3. On/off button               |
| 4. Mode selection button | 5. Power Mode indicator                    | 6. Everyday mode indicator     |
| 7. Gentle mode indicator | 8. Pulse mode indicator                    | 9. Baby mode indicator         |
| 10. Charging icon        | 11. Water container                        | 12. Rotating nozzle attachment |
| 13. Charging socket      | 14. Water filler opening for the container |                                |

Description of washing nozzles:

- A. Nozzle for daily use. For daily rinsing of interdental spaces and gum lines.
- B. Orthodontic nozzle. Nozzles for those with braces, trays or other fixed items in the mouth.
- C. Periodontal nozzle. Thin, soft tips for gently rinsing gum pockets, used, for example, for periodontal disease or after surgical procedures.
- D. Plaque nozzle. Tips with bristles or a special tip for removing plaque in hard-to-reach areas (such as the backs of molars).
- E. Tongue cleaning nozzle. For cleaning the tongue of bacterial plaque and food debris, which helps reduce bad breath, among other things.

Before first use:

- 1. Remove all packaging and marketing materials.
- 2. Using a damp clean cloth, clean all surfaces of the device.
- 3. Remove the water container after by turning the container (11) counterclockwise and rinse.
- 4. Put it back in and turn it to the right to secure the container.

Using the device:

- 1. Insert the nozzle of your choice into the rotating nozzle mount (12). You can use the swivel mount to position the nozzle at the desired angle.
- 2. Open the water filler opening and pour in clean water.
- 3. Close the filler.
- 4. Use the mode selection button (4) to select the mode that suits you.
- 5. Lean over the sink or bowl as water will flow from your mouth and the device during use.
- 6. Press the On/Off button (3). Note the water pressure may be high when turned on you can adjust the mode using the mode selection button (4).
- 7. Point the nozzle of the device at your teeth and interdental spaces. Do not point the nozzle directly at the gums as this may cause irritation and bleeding. The recommended angle of use is 90 degrees to the gum line. See illustrations page 4.
- 8. Slowly and methodically move the nozzle along the tooth line and between teeth to rinse away food debris and

accumulated bacterial plaque. It is recommended to clean the front and back of the upper jaw first and then the inside of the jaw and then the outside.

Note: Use only clean water in the water container. Mouthwash is not recommended, as it can irritate injured gums.

NOTE: With occasional use of the irrigator, the gums may bleed, more frequent use will reduce gum bleeding.

Charging the device:

1. Pull back the charging port flap (13).
2. Insert the power cord into the power port.
3. The charging icon (10) will turn red. When fully charged, it will change to green.

NOTE: Do not use the device during charging.

4. When charging is complete, close the charging port flap.

Charging time is 3-6 hours depending on the voltage of the charger.

Cleaning and maintenance:

1. Pour the water out of the water container after each use.

2. After each use, clean all nozzles used that day.

Technical specifications:

Charging power supply: 5V DC 1A

Battery capacity: 1800 mAh

Battery Type: Lithium-ion

Spare nozzle kit model: AD 2185.1

Water container capacity: 350 ml

**For the sake of the environment. Information about waste electrical and electronic equipment** In accordance with Article 13(1) and (2) of the Law of September 11, 2015 on waste electrical and electronic equipment, we inform you about the proper handling of waste electrical and electronic equipment:

1. It is forbidden to place waste electrical and electronic equipment together with other waste - this is confirmed by marking in the form of a "crossed out basket", ordering selective collection of this type of waste.

2. Electrical and electronic equipment may contain hazardous substances, mixtures and components that, when released into the environment, can cause serious danger to the health and life of humans and living organisms. They can lead to numerous health ailments, such as vision, hearing and speech disorders, and can also damage the kidneys, liver and heart, and cause skin diseases. Harmful substances, can also have adverse effects on the respiratory and reproductive systems and lead to cancerous changes. Consumption of plants growing on contaminated soils, and products made from them, can risk the aforementioned health effects.

3. Waste electrical and electronic equipment should be handed over only to authorized collection points, a list of which should be included on each Municipality's website.

4. The household plays an important role, in contributing to the reuse and recovery, including recycling, of waste equipment. It also plays a key role in the waste electrical and electronic equipment management system due to the possibility of direct delivery to authorized collection points and the elimination of undesirable social habits resulting in leaving waste equipment in places not intended for it.

**In addition**, Donate waste electrical and electronic equipment at the point of delivery. A distributor, when supplying equipment intended for households to a purchaser, is obliged to take back used household equipment free of charge at the place of delivery of that equipment, provided the used equipment is of the same type and performed the same functions as the supplied equipment.

Cardboard packaging and polyethylene (PE) bags should be placed in the appropriate containers for separate collection of municipal waste according to their description. If there are batteries in the equipment, remove them and separately dispose of them at a collection and storage facility.

**Do not dispose of the device in the municipal waste container!!!**

**Service** If you wish to purchase spare parts or report any complaints, please contact directly with the seller who issued the receipt.

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